

IMPORTANT: READ ENTIRE OPERATING MANUAL BEFORE USE.


1.



Attach the coiled twin tubing to the housing.

Screw the metal fitting on the tube into the metal fitting on the main housing. Push the white tube with no fitting into the plastic fitting on the housing. Insure each of the fittings is screwed on snugly.


2.

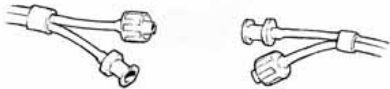

Connect the hand piece to the main housing.

Insure each connector is twisted snugly.


3.

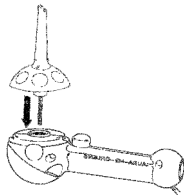

Push individual tip on.


Tip can be rotated to direct the spray.

4.

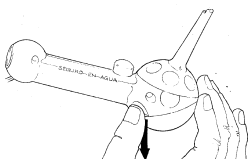

Snap cup collar off.

5.

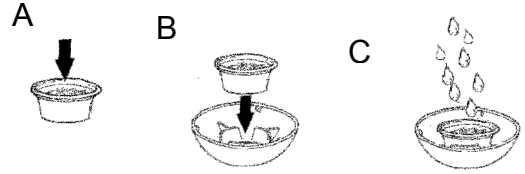


A. Fill cleaner cup with 1 tsp. cleaner
More for drier and more aggressive cleaning.

B. Place cup into collar.


C. Fill cup with water or hydrogen peroxide.
Hydrogen peroxide will give a better cleaning while also bleaching teeth.

6.

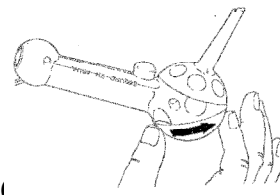

Holding cup & collar in one hand and hand piece in other, firmly snap hand piece on collar. Twist collar to set seal.

or

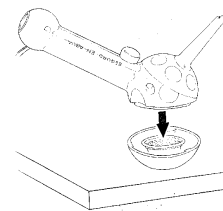

Place cup & collar on the counter. Firmly snap hand piece on collar. Twist collar to set seal.

7.



Shake hand piece until dental cleaner becomes a milkshake consistency (shake periodically throughout the cleaning process).

- Cleaner will not entirely dissolve.
- It is normal for a slight amount of cleaner to remain in the cup, particularly when the hand piece is not shaken during the cleaning process.
- The more cleaner and less water that is initially put in the cup will also result in slightly more cleaner being left in the cup. Additional liquid can be added when the liquid runs out to extend the cleaning.


8.



Press the POWER switch on the main housing.

This switch can be left on if desired.

9.

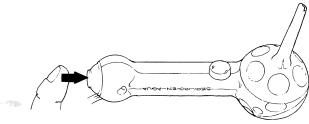

Press button on the bottom of hand piece to turn appliance ON.

Air ONLY will emit from the tip.

10.




Place mouth over sink.
(or hand piece can be taken into shower.)

Place tip in mouth.

Push or tap button under the tip to emit cleaner.

11.

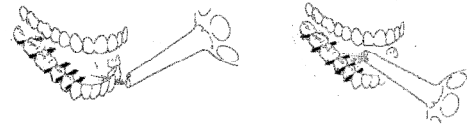

Clean teeth.


Cleaner is dispensed with air each time the button under the tip is pushed.

A light tap will emit a small amount of cleaner. A harder or longer push will emit a larger amount of cleaner.

Do not keep the button under the tip depressed continuously or it will emit all the cleaner at once.


It is best to start at the back top teeth, working around as if painting each tooth and end in the front bottom tongue-side using the slurry already in the mouth. Spend at least 2 minutes cleaning teeth, longer if periodontal disease exists.

12.

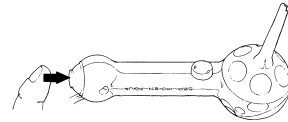

Press bottom button to shut OFF.

13.

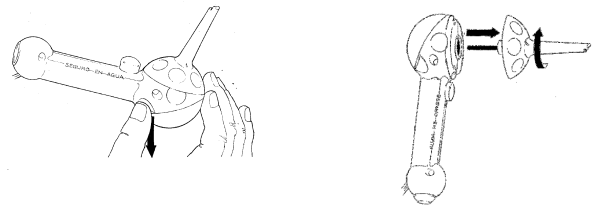

Remove cup and rinse.

Rinse entire hand piece under running water. Press button near tip while rinsing the back to clear lines.

Grasp tip, twist slightly to break seal, pull from hand piece and rinse.